

Multi RS-60, Programmable rotator




DESCRIPTION

Programmable Rotator **Multi RS-60** performs several motion types in one module. Microprocessor control allows performing not only Vertical rotation (1) of the platform, but also Reciprocal rotation (2) as well as Vibration (3). These three motion types can be performed separately, pairwise and in cycles, periodically repeating the sequence of three motion types. Multi-rotation option of Biosan instruments substantially expands possibilities and enhances efficiency of sample preparation for the examined materials and allows adjusting the mixing procedure according to the individual tasks.

Programmable Rotator can be used for variety of applications in modern life science laboratories: for hybridization reactions, cell growing, soft extraction and homogenisation of biological components in solutions, as well as for reactions of binding and washing of magnetic particles.

During the stop mode the platform does not perform extra turnover and stops in a horizontal plane. Additionally there is a possibility of setting the plane with respect to which platform oscillates vertically or horizontally.

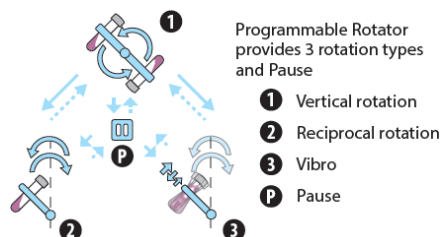

CAT. NUMBER

| Including PRS-48 platform | Including PRS-48 platform |
|---------------------------|---|
| BS-010118-AAI | 230VAC 50/60Hz Euro plug |
| BS-010118-AAQ | 230VAC 50/60Hz UK plug |
| BS-010118-AA4 | 230VAC 50/60Hz AU plug |
| BS-010118-AA2 | 100VAC 50/60Hz US plug, 120VAC 60Hz US plug |
| BS-010118-DK | IQ OQ document |
| BS-010118-EK | PQ document |

SPECIFICATIONS

| | |
|---|--|
| Digital time setting | 1 min–24 hrs / non–stop (increment 1 min) |
| Timer sound signal | + |
| Vertical rotation speed control range | 1–100 rpm (increment 1 rpm) |
| Vertical rotation movement | overhead, 360° |
| Vertical rotation time setting range | 0 - 250 sec |
| Reciprocal rotation speed control range | 1–100 rpm (increment 1 rpm) |
| Reciprocal rotation tilt angle range | 1° - 90° (increment 1°) |
| Reciprocal rotation time setting range | 0 - 250 sec |
| Vibro rotation tilt angle range | 0° - 5° (increment 1°) |
| Vibro rotation Pause / Vibro time setting range | 0 - 5 sec |
| Maximum load | 0.8 kg |
| Overall dimensions (W×D×H) | 430×230×230 mm |
| Weight | 3.8 kg |
| Input current/power consumption | 24 V, 750 mA / 18 W |
| External power supply | Input AC 100–240 V; 50/60 Hz; Output DC 24 V |

ACCESSORIES



PRS-14
BS-010118-BK
platform

14 tubes 20-30 mm diameter
(50 ml tubes).



PRS-48
BS-010118-CK
platform

48 tubes 10-16 mm diameter
(1.5 ml-15 ml tubes).



PRS-8/22
BS-010118-AK
platform

8 tubes 20-30 mm diameter (50
ml tubes), 22 tubes 10-16 mm
diameter (1.5 ml-15 ml tubes).